


Balanced Aire Non-Powered Self Adjusting Convertible Mattress



HCPCS: E0373



FIGURE A

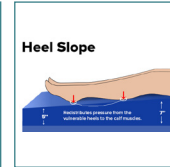

INSET B



INSET C

- Versatile reactive mattress replacement system with self adjusting air/foam cells provides outstanding pressure redistribution and shear/friction reduction
- Designed to prevent, treat and heal pressure ulcers in the home or long term care setting
- Add an alternating pressure pump and the Balanced Aire easily converts to a dynamic powered system
- 2.5" Hybrid Visco foam top layer provides immersion, premium comfort and heat dissipation
- High density foam perimeter provides firm edge support and facilitates safe patient transfers
- Heel slope redistributes pressure away from vulnerable heels
- Self-adjustment can be open cell or closed cell to meet individual needs
- Self-Adjusting air cells recognize body weight and pressure and adjust to reduce interface pressure and improve patient comfort
- Eight non-circular cylinders prevent "rolling" or "gapping"
- Removable 4 way stretch PU cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- Arrives patient ready; easy to deliver, set up and transport
- Foam and cover are antimicrobial, antifungal, and latex free
- Meets flammability standards CFR1632 and CRF1633

PRODUCTS

Description

BA9600-NP	Non-Powered Self Adjusting Air/Foam Mattress, E0373, 1/ea
BA9600-P	Powered alternating Pressure Air/Foam Mattress, E0277, 1/ea
BA9600-NP-84	Non-Powered Self Adjusting Air/Foam Mattress, E0373, 1/ea
BA9600-P-84	Powered Alternating Pressure Air/Foam Mattress, E0277, 1/ea
BA9600-NP-42	Non-Powered Self Adjusting Air/Foam Mattress, E0373, 1/ea
BA9600-P-42	Powered Alternating Pressure Air/Foam Mattress, E0277, 1/ea

SPECIFICATIONS

Description

Control Unit Airflow	When Used w/Power Source (P) 4 LPM
Control Unit Cycle Time	When Used w/Power Source (P) 5 Minutes
Control Unit Type	When Used w/Power Source (P) Compressor
Control Unit Alarms	N/A
Mattress Cover Material	Fluid-Resistant Stretch
Weight capacity	35" W- 500 lbs., 42" W- 700 lbs.
Warranty	5 Year Limited Mattress 1 Year Control Unit

ACCESSORIES

Description

14005E	Variable Pressure Pump Only for 14001E (Replaces Mason # 2002-ADJ), 1/ea
BA9600-P-42-TC	Replacement Cover for 42", 1/ea
BA9600-C-NP3584	Replacement cover for 35" (W) x 84" (L), 1/ea
BA9600-C-NP3580	Replacement cover for 35" (W) x 80" (L), 1/ea